



Alumina Disc
Diameter: 14mm



Sapphire Window
Diameter: 9.5mm
Thickness: 0.5mm



Graphite Crucible Body top view

Inner diameter: 14.55mm

Outer diameter: 20.30mm

Height: 6.5mm

Bottom thickness: 2mm





Crucible Body Bottom view
Diameter: 20.30mm





Graphite Lid Oblique View

Lower D: 14.55mm

Upper D: 13.70mm

Height 1: 4mm

Height 2: 6mm

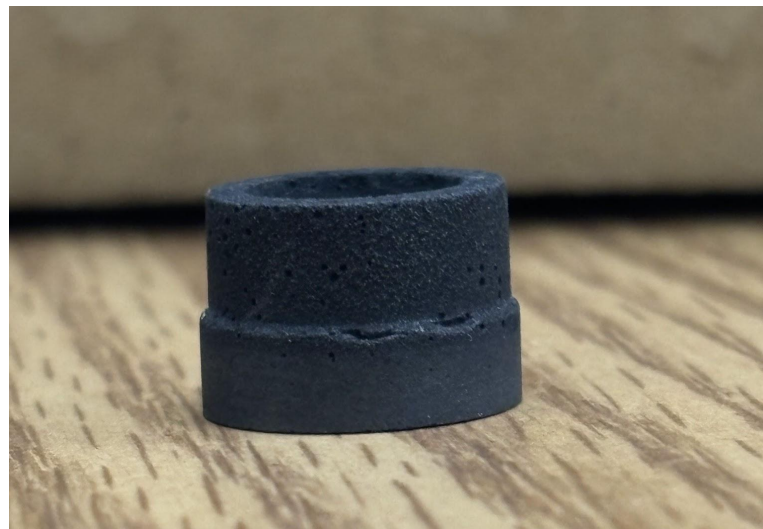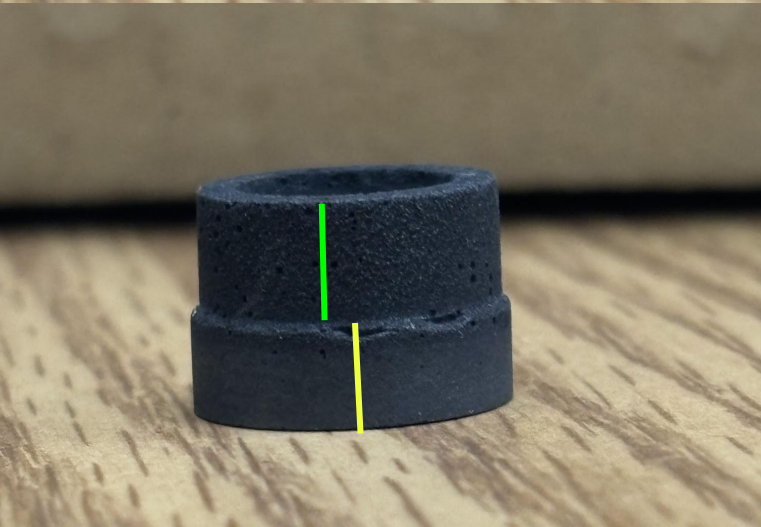




Graphite Lid Top View

Outermost Diameter: 13.70mm

2nd Diameter: 9.75mm

3rd Diameter: 5mm

4th Diameter: 3.5mm

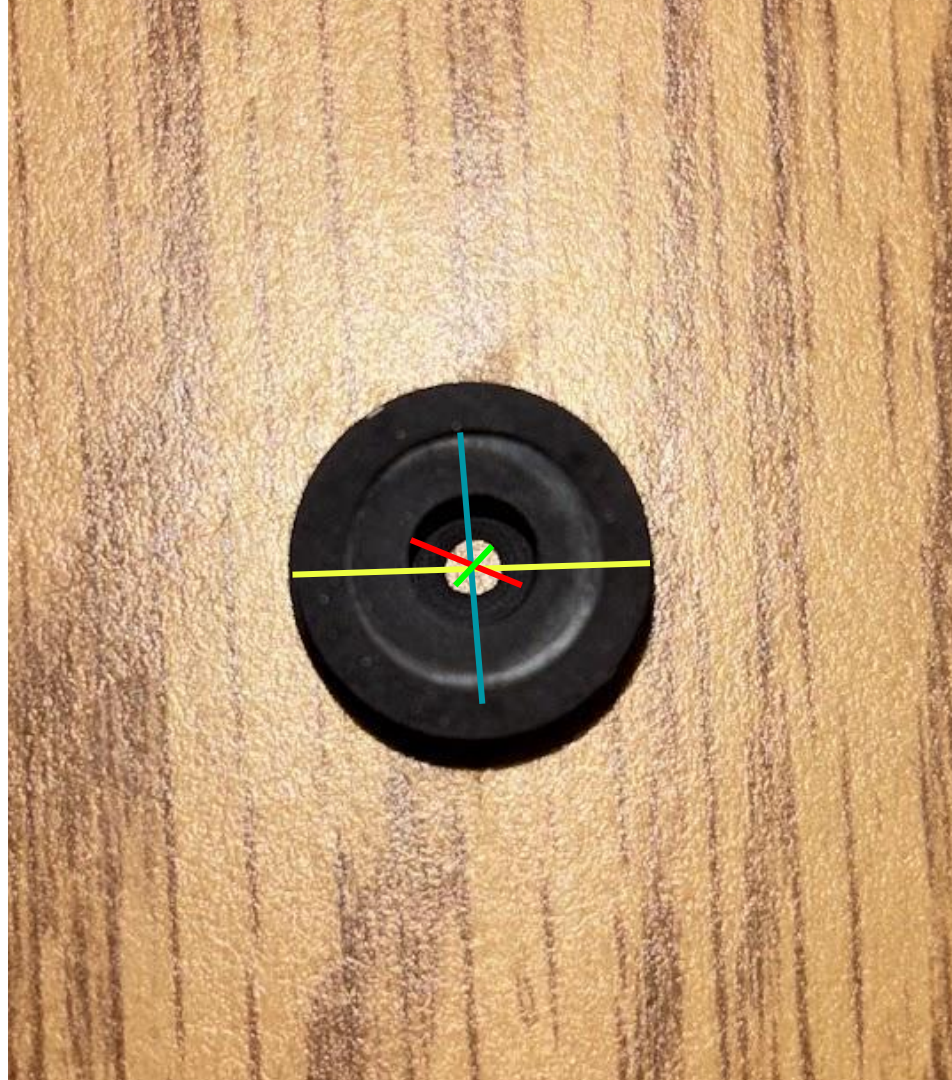


Graphite Lid Bottom View

Inner Diameter: 3.5mm

Outer Diameter: 14.55mm



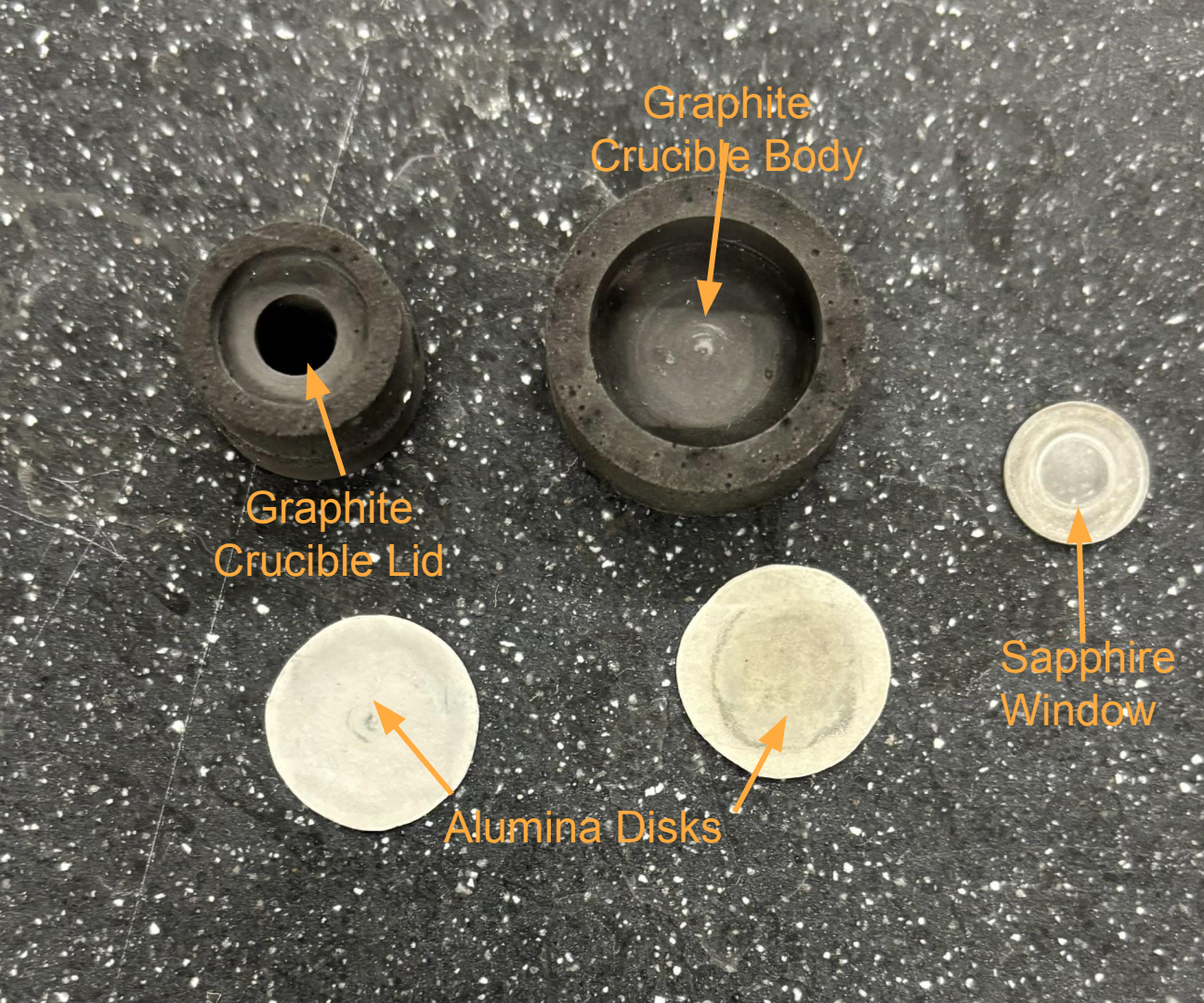


Graphite
Crucible Body

Graphite
Crucible Lid


Alumina Disks

Sapphire
Window

Caption: All parts of
the crucible layed
next to each other.

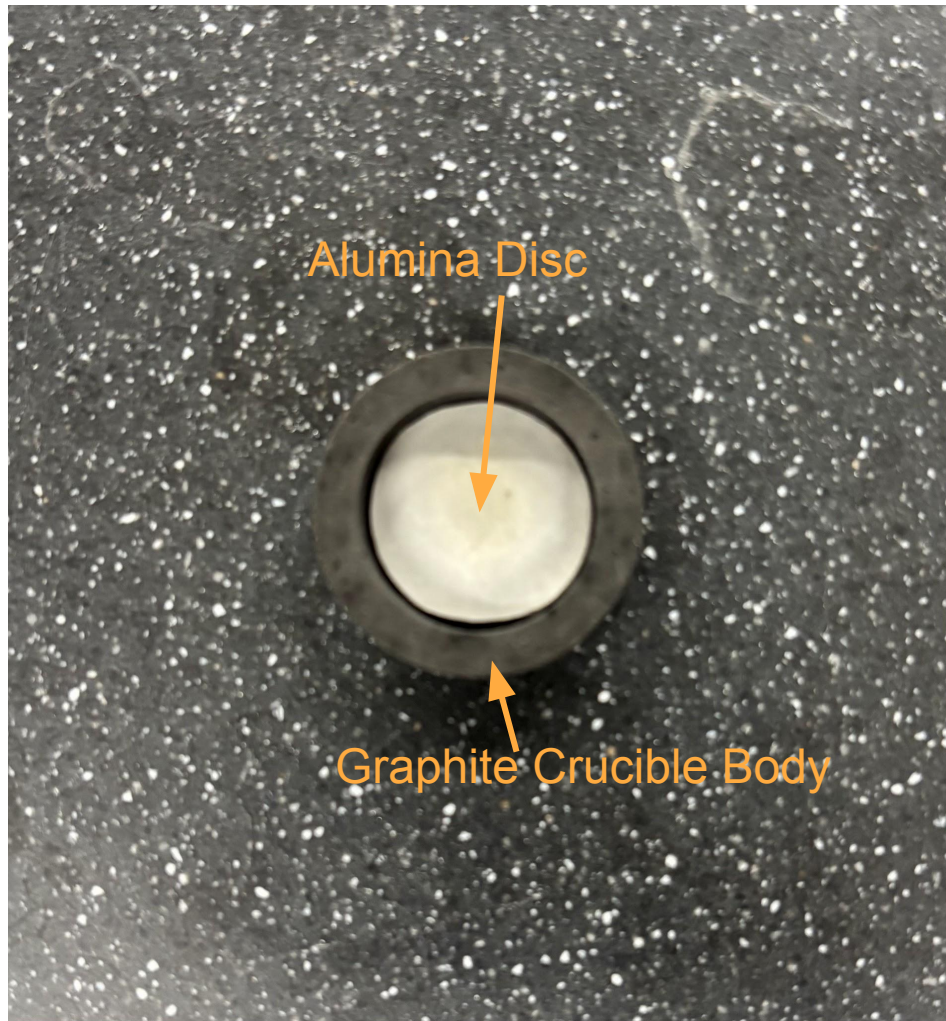

Caption: An alumina sheet rests in the bottom of the crucible. The sample is placed on top of the alumina sheet, then covered with the second alumina sheet.

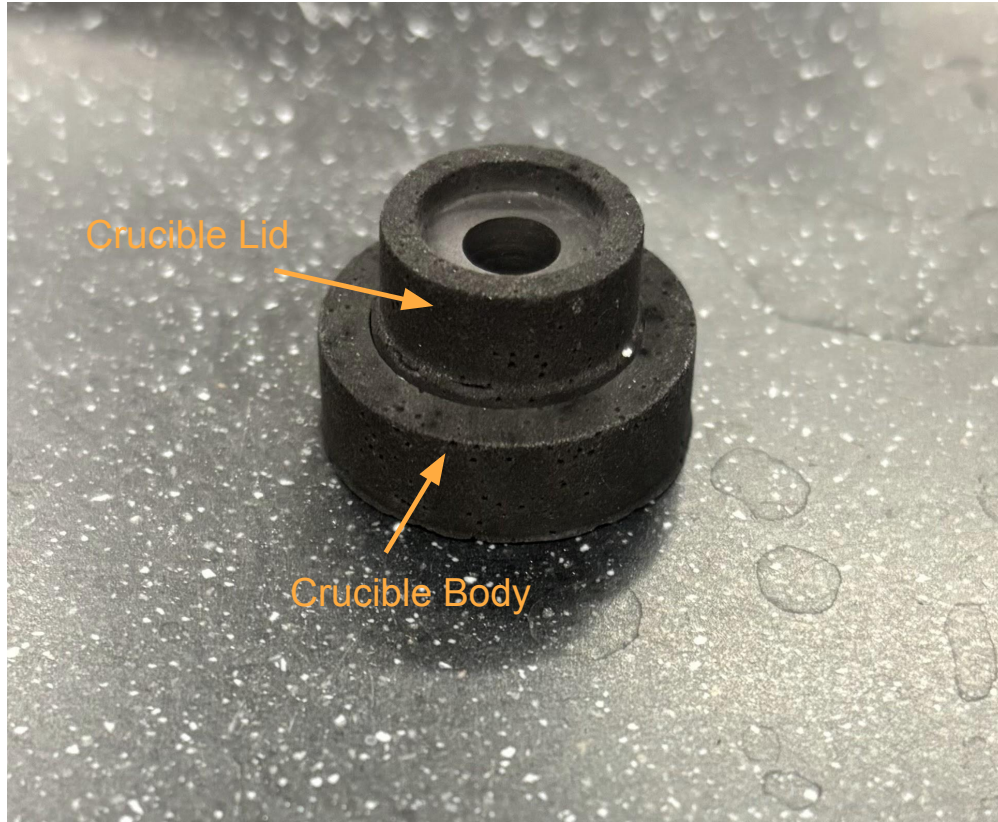

Caption: The lid is placed on top of the alumina-sample stack



Caption: The sapphire window is placed on top of the lid, allowing the pyrometer to read the sample temperature.